



THE PRINCIPLE OF THE PRESSOTHERAPY



The pressotherapy, a process of drainage and reshaping, uses a set of accessories such as boots, cuffs or belt to relieve the targeted part of the body.


These accessories are endowed with compartments which fill with air in different rates, creating infinity of small pressures on the treated zone. The pressotherapy stimulates the blood circulation. Now, It is used to remodel the silhouette and for its draining effect on the cellulite.

PNEUMATIC LYMPHATIC DRAINAGE: PRESSOFAST TECHNOLOGY

The benefactions of the lymphatic drainage are not to be demonstrated any more. PRESSOFAST lavishes soft mechanical palpitations. For example the boots, this one are divided in compartments. The first compartment, the one who surrounds the foot, is inflated by a pressure induced by PRESSOFAST. Once the pressure is

established in the second compartment, the pressure exercised in the initial compartment increases gradually. The difference of pressure so created causes the migration of the lymph in the high-pressure compartment to the low-pressure compartment. In a similar way, the following compartments

are put under pressure. During the treatment, the modification of the pressure depends on the volume of the treated member. Consequently, PRESSOFAST insures efficiency throughout the treatment, do not damage the superficial vessels thanks to its technique of soft palpations, and is the ideal support for the reshaping.

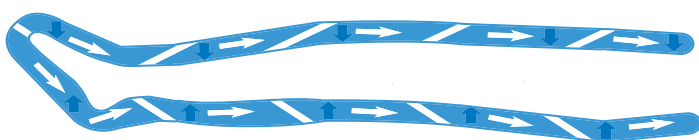

MIGRATION OF THE LYMPH IN THE FOOT TO THE DIRECTION WANTED THANKS TO THE PROGRESSIVE PRESSURE EXERCISED ON THE COMPARTMENTS OF THE BOOT



Copyright ©DERMEO* - Illustrations and information non-contractual. The products and treatments offered by DERMEO* or its distributors may be subject to specific regulations.

CHARACTERISTICS OF PRESSOFAST

PRESSOFAST is a pneumatic lymphatic drainage device with a pressure automatic modulation, fully configurable. This is the easiest device to use thanks to its easiness of its regulation (treatment time per room, treatment time, pressure ...).

Dimension in mm (H x W x D)	L400mm x l310mm x H110mm
Weight	6.5 kg
Range of pressure	20 - 200 mB (15-150mmHg)
Consumption	44 w
Intensity/Frequency	220V/50 V
Sound intensity	49dB < I < 52dB
Easy of use	Thanks to its forward technology, PRESSOFAST is very noiseless

Mode of the pressure regulation	Double Gradient of Pressure Automatic
Adjustable durations	Treatment time - Deflating - Rest and processing time
Inflating of elements	2 rows of 5 exits which allow the inflating of 1 or 2 elements simultaneously
Composition of elements	Outside : Nylon Inside : Polyurethane Closure : Thermo-welded elements

ACCESSORIES



CUFF MODEL WITH 5 COMPARTMENTS 58 X 77 CM



ABDOMINAL MODEL WITH 5 COMPARTMENTS 56 X 37 CM



BOOT MODEL WITH 5 COMPARTMENTS 80 X 92 CM

PRESSOFAST APPLICATIONS

PRESSOFAST applications are numerous in particular to relieve the heavy legs and reduce the cellulite.

HEAVY LEGS

The pressotherapy targets the heavy legs, indeed, it strengthens the blood return which materializes by a reduction of the inflation and the decrease of the sensation of heavy legs. A feeling of well-being and comfort is felt from the first session.

CELLULITE

The device PRESSOFAST enables to reduce the orange skin effect because it leads a decrease of water retention. To perpetuate the obtained results, the sessions can be accompanied with a light diet.



Dermeo
Paris

©DERMEO® - Illustrations and information non-contractual. The products and treatments offered by DERMEO® or its distributors may be subject to specific regulations. Also it will suit the user or purchaser of the equipment to check with its distributor and competent authorities on the regulations applicable in the country of use as well as the specific conditions of use of such equipment.